# **Anti-CGAS antibody**

Catalog Number: 176598



#### **Product name**

Anti-CGAS antibody

## **Specificity**

Human, Mouse, Rat

## **Antibody description**

Rabbit polyclonal antibody to CGAS

## Preparation

This antigen of this antibody was synthetic peptide within c-terminal human cgas.

#### **Formulation**

Liquid, 1\*TBS (pH7.4), 0.2% BSA, 50% Glycerol. Preservative: 0.05% Sodium Azide.

## **Storage**

Store at +4°C after thawing. Aliquot store at -20°C. Avoid repeated freeze / thaw cycles.

### Clonality

Polyclonal

#### Ig Type

Rabbit IgG

#### **Applications**

WB, ICC, IHC-P

### **Dilutions**

WB: 1:500-1:1,000

ICC: 1:50-1:100

IHC-P: 1:100-1:500

## Validations


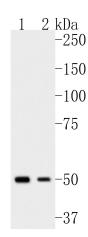


Fig1:; Western blot analysis of cGAS on different lysates. Proteins were transferred to a PVDF membrane and blocked with 5% BSA in PBS for 1 hour at room temperature. The primary antibody ( 1/500) was used in 5% BSA at room temperature for 2 hours. Goat Anti-Rabbit IgG - HRP Secondary Antibody (HA1001) at 1:5,000 dilution was used for 1 hour at room temperature.; Positive control:; Lane 1: MCF-7 cell lysate; Lane 2: Mouse lung

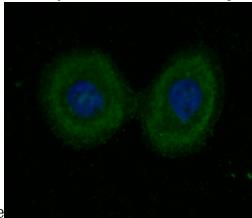

tissue lysate

Fig2:; ICC staining of cGAS in A431 cells (green). Formalin fixed cells were permeabilized with 0.1% Triton X-100 in TBS for 10 minutes at room temperature and blocked with 1% Blocker BSA for 15 minutes at room temperature. Cells were probed with the primary antibody (1/100) for 1 hour at room temperature, washed with PBS. Alexa Fluor®488 Goat anti-Rabbit IgG was used as the secondary antibody at 1/1,000 dilution. The nuclear counter stain is DAPI (blue).

# **Anti-CGAS antibody**

Catalog Number: 176598



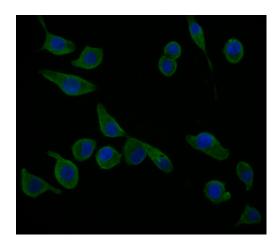


Fig3:; ICC staining of cGAS in LOVO cells (green). Formalin fixed cells were permeabilized with 0.1% Triton X-100 in TBS for 10 minutes at room temperature and blocked with 1% Blocker BSA for 15 minutes at room temperature. Cells were probed with the primary antibody (1/100) for 1 hour at room temperature, washed with PBS. Alexa Fluor®488 Goat anti-Rabbit IgG was used as the secondary antibody at 1/1,000 dilution. The nuclear counter stain is DAPI (blue).

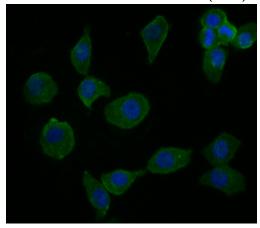


Fig4:; ICC staining of cGAS in SKOV-3 cells (green). Formalin fixed cells were permeabilized with 0.1% Triton X-100 in TBS for 10 minutes at room temperature and blocked with 1% Blocker BSA for 15 minutes at room temperature. Cells were probed with the primary antibody (1/100) for 1 hour at room temperature, washed with PBS. Alexa Fluor®488 Goat anti-Rabbit IgG was used as the secondary antibody at 1/1,000 dilution. The nuclear counter stain is DAPI (blue).

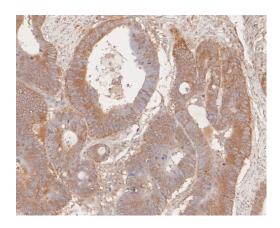


Fig5:; Immunohistochemical analysis of paraffinembedded human colon carcinoma tissue using anti-cGAS antibody. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0) for 20 minutes. The tissues were blocked in 5% BSA for 30 minutes at room temperature, washed with ddH; 2; O and PBS, and then probed with the primary antibody (1/200) for 30 minutes at room temperature. The detection was performed using an HRP conjugated compact polymer system. DAB was used as the chromogen. Tissues were counterstained with hematoxylin and mounted with DPX.

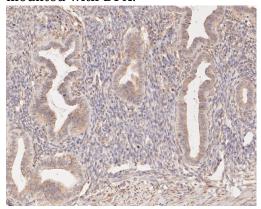


Fig6:; Immunohistochemical analysis of paraffinembedded human uterus tissue using anti-cGAS antibody. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0) for 20 minutes. The tissues were blocked in 5% BSA for 30 minutes at room temperature, washed with ddH; 2; O and PBS, and then probed with the primary antibody (1/200) for 30 minutes at room temperature. The detection was performed using an HRP conjugated compact polymer system. DAB was used as the chromogen.

# **Anti-CGAS** antibody

Catalog Number: 176598



Tissues were counterstained with hematoxylin and mounted with DPX.

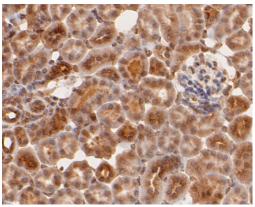


Fig7:; Immunohistochemical analysis of paraffinembedded mouse kidney tissue using anti-cGAS antibody. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0) for 20 minutes. The tissues were blocked in 5% BSA for 30 minutes at room temperature, washed with ddH; 2; O and PBS, and then probed with the primary antibody (1/200) for 30 minutes at room temperature. The detection was performed using an HRP conjugated compact polymer system. DAB was used as the chromogen. Tissues were counterstained with hematoxylin and

mounted with DPX.

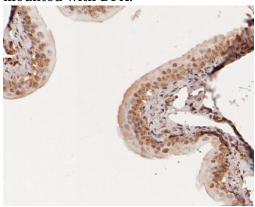


Fig8:; Immunohistochemical analysis of paraffinembedded rat bladder tissue using anti-cGAS antibody. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0) for 20 minutes. The tissues were blocked in 5% BSA for 30 minutes at room temperature, washed with ddH; 2; O and PBS, and then probed with the primary antibody (1/50) for 30 minutes at room temperature. The detection was performed using an HRP conjugated compact polymer system. DAB was used as the chromogen. Tissues were counterstained with hematoxylin and mounted with DPX.