# **Anti-DYSF antibody**

Catalog Number: 176542



#### **Product name**

Anti-DYSF antibody

# **Specificity**

Human, Mouse

# **Antibody description**

Rabbit monoclonal antibody to DYSF

# Preparation

This antigen of this antibody was synthetic peptide within human dysferlin aa 110-150.

#### **Formulation**

Liquid, 1\*TBS (pH7.4), 0.05% BSA, 40% Glycerol. Preservative: 0.05% Sodium Azide.

# **Storage**

Store at +4°C after thawing. Aliquot store at -20°C or -80°C. Avoid repeated freeze / thaw cycles.

#### **Clonality**

Monoclonal

#### Ig Type

Rabbit IgG

#### **Applications**

WB, ICC/IF, IHC-P

### **Dilutions**

WB: 1:1,000-1:5,000

ICC/IF: 1:50-1:200

IHC-P: 1:50-1:200

#### **Validations**

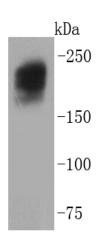


Fig1:; Western blot analysis of Dysferlin on skeletal muscle tissue lysates. Proteins were transferred to a PVDF membrane and blocked with 5% BSA in PBS for 1 hour at room temperature. The primary antibody ( 1/500) was used in 5% BSA at room temperature for 2 hours. Goat Anti-Rabbit IgG - HRP Secondary Antibody (HA1001) at 1:40,000 dilution was used for 1 hour at room temperature.

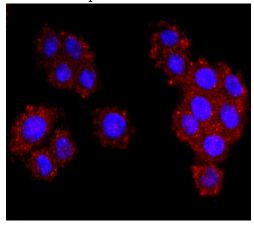


Fig2:; ICC staining of Dysferlin in SW480 cells (red). Formalin fixed cells were permeabilized with 0.1% Triton X-100 in TBS for 10 minutes at room temperature and blocked with 1% Blocker BSA for 15 minutes at room temperature. Cells were probed with the primary antibody (1/50) for 1 hour at room temperature, washed with PBS. Alexa Fluor®594 Goat anti-Rabbit IgG was used as the secondary antibody at 1/1,000 dilution. The nuclear counter stain is DAPI (blue).

# **Anti-DYSF antibody**

Catalog Number: 176542



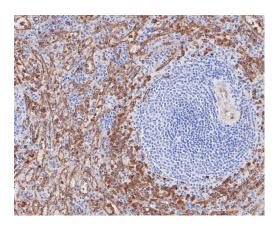


Fig3:; Immunohistochemical analysis of paraffinembedded human spleen tissue using anti-Dysferlin antibody. The section was pre-treated using heat mediated antigen retrieval with Tris-EDTA buffer (pH 8.0-8.4) for 20 minutes. The tissues were blocked in 5% BSA for 30 minutes at room temperature, washed with ddH; 2; O and PBS, and then probed with the primary antibody (1/50) for 30 minutes at room temperature. The detection was performed using an HRP conjugated compact polymer system. DAB was used as the chromogen. Tissues were counterstained with hematoxylin and mounted with DPX.

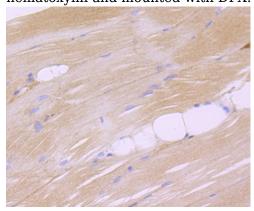


Fig4:; Immunohistochemical analysis of paraffinembedded mouse skeletal muscle tissue using anti-Dysferlin antibody. The section was pretreated using heat mediated antigen retrieval with Tris-EDTA buffer (pH 8.0-8.4) for 20 minutes. The tissues were blocked in 5% BSA for 30 minutes at room temperature, washed with ddH; 2; O and PBS, and then probed with the primary antibody (1/50) for 30 minutes at room temperature. The detection was performed using an HRP conjugated compact polymer system. DAB was used as the

chromogen. Tissues were counterstained with hematoxylin and mounted with DPX.

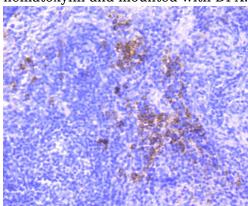


Fig5:; Immunohistochemical analysis of paraffinembedded mouse spleen tissue using anti-Dysferlin antibody. The section was pre-treated using heat mediated antigen retrieval with Tris-EDTA buffer (pH 8.0-8.4) for 20 minutes. The tissues were blocked in 5% BSA for 30 minutes at room temperature, washed with ddH; 2; O and PBS, and then probed with the primary antibody (1/50) for 30 minutes at room temperature. The detection was performed using an HRP conjugated compact polymer system. DAB was used as the chromogen. Tissues were counterstained with hematoxylin and mounted with DPX.

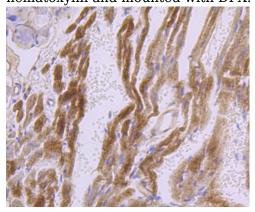


Fig6:; Immunohistochemical analysis of paraffinembedded mouse heart tissue using anti-Dysferlin antibody. The section was pre-treated using heat mediated antigen retrieval with Tris-EDTA buffer (pH 8.0-8.4) for 20 minutes. The tissues were blocked in 5% BSA for 30 minutes at room temperature, washed with ddH; 2; O and PBS, and then probed with the primary antibody (1/50) for 30 minutes at room temperature. The detection

# **Anti-DYSF antibody**

Catalog Number: 176542



was performed using an HRP conjugated compact polymer system. DAB was used as the chromogen.

Tissues were counterstained with hematoxylin and mounted with DPX.